

DONECO

King Side Heating Pad, Electric Foot Warmer with 4 Temperature Settings and Fast-Heating Technology



DN0508

Instruction Manual


PRODUCT SPECIFICATIONS

Pad Size: 22" x 22"
Power Supply: AC 120V, 60Hz, 55W
Pad Material: Micromink
Auto Shut-Off Time: 120 minutes

4 HEAT SETTINGS:

Warm-green light, temperature range (~110°F)
Low-yellow light, temperature range (~120°F)
Medium-orange light, temperature range (~130°F)
High-red light, temperature range (~140°F)

PRODUCT FEATURES

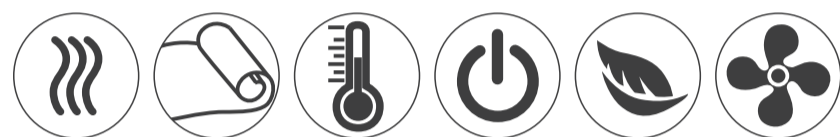


- 1 Fast heating:** Heats in seconds to provide fast pain relief to aching muscles.
- 2 The extra large pad is sized at 22" x 22":** Perfect for relief on many different muscle groups.
- 3 4 Temperature Settings:** Digital LED controller adjusts heat easily with 4 temperature settings to accommodate your comfort level.
- 4 Microplush Material:** Thick microplush fibers offer maximum comfort for daily use.
- 5 2-Hour Auto Shut-Off:** Conserves energy and prevents burns caused by excessive heating.
- 6 Machine Washable:** Helps maintain beauty and softness for years to come.

2

USING DN0508

Attach the Pad Connector to the Control Connector. Push both parts together securely to ensure a tight fit.

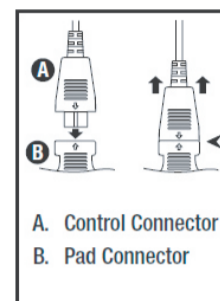

▲ NOTE: The Control Connector must be connected to the Pad Connector before plugging the unit into an electrical outlet.

- After connecting the Pad Connector the Control Connector, plug the unit into a 120V AC electrical outlet.
- To turn the DN0508 on, press the power button to select the desired heat setting:
Warm-green light, temperature range (~110°F)
Low-yellow light, temperature range (~120°F)
Medium-orange light, temperature range (~130°F)
High-red light, temperature range (~140°F)
- To turn the DN0508 off, press the power button one more time.

▲ NOTE: If the unit does not turn on when you press the power button, make sure the Control Connector is securely plugged into the Pad Connector.

USING MOIST HEAT THERAPY

In some cases, you may find that moist heat therapy provides more comfort and tension relief. When using moist heat therapy, you should protect furniture and clothing from becoming damp.

3

- Dampen one side of the DN0508 with water by patting the surface with a wet sponge or cloth or use a spray bottle to apply a fine mist. Do not thoroughly soak the DN0508.
- Do not spray or dampen the connector.



DANGER: Do not immerse any part of the DN0508 in water or place under running water when plugged in. Doing so can damage the DN0508 and cause serious personal injury.



▲ NOTE: The Pad Connector and Control Connector must be completely dry before establishing an electrical connection.

AUTO-OFF FEATURE

To conserve energy and prevent skin burns caused by excessive use, the DN0508 includes an auto-off feature. The unit will turn off automatically after two hours of continuous use. You can turn this feature off by pressing the Auto Off button until the Stay On indicator light turns red. The unit will not turn off automatically after two hours of continuous use when the Stay On light is red.

▲ NOTE: After a power outage, the DN0508 automatically resets the auto-off feature.

WASHING INSTRUCTIONS

Do Not Dry Clean! Do Not Use Bleach! Do Not Use Wringer!

- Disconnect the DN0508 from the electrical outlet.
- Detach the Pad Connector from the Control Connector.
- Machine wash the DN0508 in cold water on the gentle cycle or hand wash. Line dry or tumble dry on high heat for five minutes.
- If using an automatic dryer, remove while damp and drape to finish drying.
- Do not iron the heating pad.
- Do not plug in the DN0508 until the Pad Connector is completely dry.

4

WARRANTY

The DN0508 comes with an industry-leading 5 Year Warranty that begins on the date of purchase.

The warranty applies to the heating pad and necessary parts and labor relating thereto. The warranty does not apply to damage resulting from failure to follow the operating instructions, accidents, abuse, alterations, or disassembly by unauthorized individuals.

ARE YOU 100% SATISFIED?

If you have product questions or concerns, don't hesitate to contact us directly at support@myfrieq.com

DONECO has an award-winning warranty, exchange, and customer service program that guarantees hassle-free solutions to any issue you might have within 24 hours!

5

DONECO

Tel: 1-909-595-1200 (US)
(Mon.9:00-17:00 PST)